


# MODEL SU105UE

Compact Ultra-Low Temperature Freezer

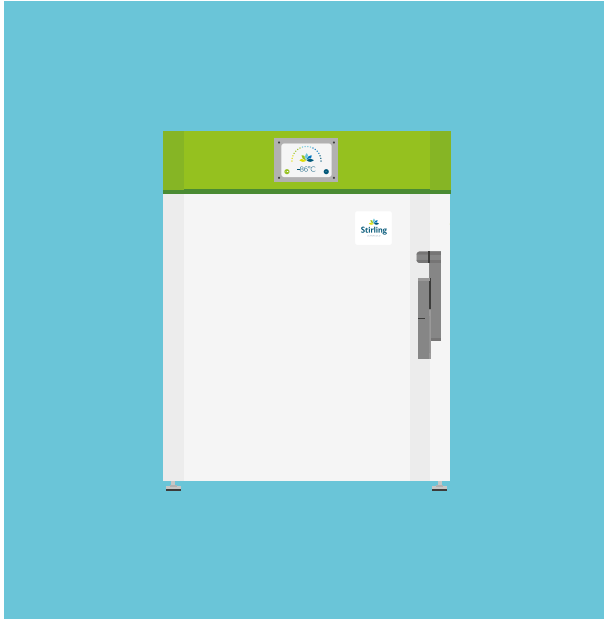

## Surprising ULT storage capacity in a space-saving design

- Advanced Stirling engine reliability maximizes sample storage in a compact footprint without sacrificing performance
- Upright design can be stacked or placed undercounter for convenient access when space is limited
- Perfect for controlled temperature transitioning within the widest range of set point temperatures, -86°C to -20°C
- Lowest energy consumption in its class for greater operational savings, <4 kWh/day at -80°C
- Universal power adaptation accepts electrical input from any outlet, around the world
- Optional SenseAnywhere wireless temperature logging allows for remote monitoring via cloud-based software
- Environmentally friendly design uses only 100% natural refrigerants for cooling



Two SU105UE freezers shown stacked. Stacking adapter kit sold separately.  
\*Stacking the SU105UE may void UL 471 certification.


## Make a big impact across your organization with the compact Stirling SU105UE

- Designed for ultra-low power consumption and maximal energy savings
- Virtually eliminates any need for maintenance with breakthrough Stirling free-piston engine
- Provides secure individual storage that is convenient to the benchtop
- Uses 100% natural refrigerants for reduced climate impact
- Industry-leading 7-year warranty for engine and thermosiphon, with two years of included parts and labor coverage\*

*\*Labor warranty coverage only available in the U.S. and Canada*

# MODEL SU105UE Specifications

## → Application, rating and electric data

<b>Application</b>	Storage of general (non-flammable) laboratory materials
<b>Storage Volume</b>	105 liters (3.7 cu.ft.)
<b>Storage Capacity</b>	72 standard 2" boxes in optional racks, available separately
<b>Temperature Range</b>	-86°C to -20°C, adjustable to 1°C increments
<b>Electric Power†</b>	100V to 240V (±10%) at either 50 or 60 Hz
<b>Power Cord Options Available</b>	Multiple options are available. Consult with your sales representative
	
<b>Maximum Power (Current)</b>	300 watts (3 amps @120V, 1.5 amps @240V) (±10%) at either 50 or 60 Hz
<b>Auto-Voltage Capability</b>	100V to 240V, 50 or 60 Hz
<b>Electric Supply Rating</b>	15 amp or greater grounded circuit
<b>Certification/Agency Listing</b>	CE, cULus
<b>Noise</b>	Advanced noise abatement, <48 dB(A) at 1 meter
<b>Indoor/Outdoor Use</b>	Indoor use only
<b>Application Environment</b>	Non-corrosive, non-flammable, non-explosive
<b>Ambient Operating Temp</b>	5°C to 35°C (41°F to 95°F)

## → Controller

<b>Interface</b>	Graphic user touchscreen interface
<b>Controller Type</b>	Microprocessor with touchscreen input and display
<b>Security</b>	Lockable door Optional PIN requirement built in
<b>Warm and Cold Alarms</b>	Fully adjustable
<b>Control Sensor</b>	RTD (PT100 Class A)
<b>Event Log</b>	All alarms, door openings
<b>Dry Contacts</b>	Normally closed, normally open, common; activated by power outage or any alarm condition (24 Vdc maximum/500mA maximum)
<b>Temperature Log</b>	30 days available graphically
<b>Battery Back-up</b>	12 hour control battery back-up for touchscreen (6V lead acid battery)
<b>Internet Connectivity</b>	Optional SenseAnywhere wireless temperature monitoring and logging

## → Refrigeration system

<b>Cooling Engine</b>	Helium charged free-piston Stirling engine with continuous modulation
<b>Heat Transport System</b>	Gravity driven thermosiphon
<b>Refrigerant</b>	R-170 (Ethane) 30–33 grams
<b>Evaporator</b>	Cold wall (inner liner)
<b>Heat Rejection</b>	Finned heat exchanger with forced air cooling Air inlet: Above freezer door, below mechanical compartment Air outlet: Underneath of rear right corner bottom
<b>Defrost Method</b>	Manual

## → Performance data

<b>Steady State Energy Use</b>	< 4 kWh/day at -80°C
<b>Pull-Down from 25°C Ambient</b>	12 hours to -80°C (Empty Cabinet)
<b>Recovery from Door Opening (ENERGY STAR® Final Test Method)</b>	40 minutes at -80°C
<b>Warm-up Profile</b>	2 hours to -60°C at -80°C (Empty Cabinet) 5.2 hours to -40°C at -80°C (Empty Cabinet)
<b>Heat Dissipation</b>	569 BTU/h (load to HVAC) at -80°C (Empty Cabinet)

## → Dimensions and construction

<b>Interior (H x D x W)</b>	533 x 432 x 457 mm / (21 x 17 x 18 in.)
<b>Exterior (H x D x W)</b>	864 x 711 x 686 mm / (34 x 28 x 27 in.)
<b>Net Weight, Empty</b>	100 kg (220 lbs.)
<b>Shipping (H x D x W)</b>	1079.5 x 838.2 x 838.2 mm / (42.5 x 33 x 33 in.)
<b>Shipping Weight</b>	121 kg (266 lbs.)
<b>Insulation</b>	High performance vacuum insulated panels and polyurethane foam using Ecomate® environmentally friendly, SNAP-compliant blowing agent
<b>Gasket Heater</b>	User programmable duty cycle
<b>Shelves</b>	1 stainless steel, fixed position (removable)

† The freezer operates on nominal 120V or 240V, 50 or 60Hz AC, over a wide voltage range from 85V to 264V. There is no need for special wiring or a 20 amp breaker on a 120V line. Heat rejection is minimal. The cooling air is drawn in over the top of the door and warm air exits at the right rear of the cabinet.



Stirling Ultracold 6000 Poston Road, Athens, Ohio 45701, USA  
T 740.274.7900 / 1.855.274.7900 | F 740.274.7901 | [www.stirlingultracold.com](http://www.stirlingultracold.com)

©2024 Stirling Ultracold, Global Cooling, Inc. All Rights Reserved.

Global Cooling technology is manufactured under U.S. and International patents. Stirling Ultracold is a trademark of Global Cooling, Inc. Specifications subject to change without notice. Refer to [www.stirlingultracold.com](http://www.stirlingultracold.com) for the latest specifications.



STIRLING POWER